

**Required Tool: (Not included)**

- 4 mm Allen wrench

**⚠ WARNING**


**⚠ CAUTION:** To avoid potential shock or other injury, turn power switch OFF and disconnect charger before conducting any assembly or maintenance procedures. Failure to follow these steps in the correct order may cause irreparable damage.

**Step 1**

Using a 4 mm Allen wrench, remove the four (4) bolts from the seat. Remove the seat from the unit.

**Step 2**

Using a 4 mm Allen wrench, remove the four (4) bolts from the seat belt bracket.

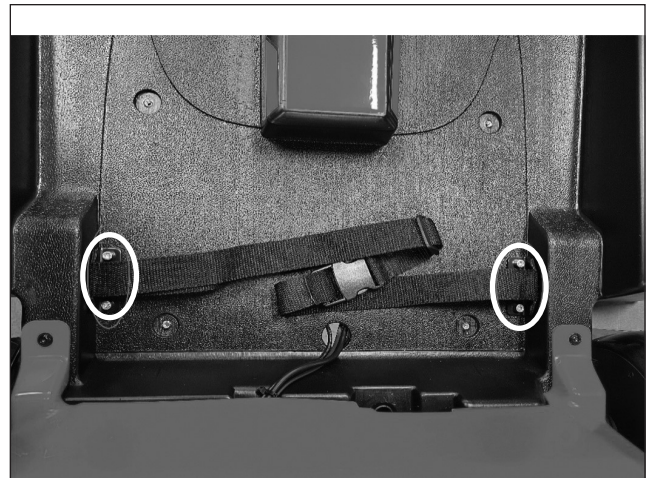

### Step 3

Remove the seat belt from the unit.



### Step 4

Reverse the steps to install the new seat belt.

**ATTENTION: Charge battery 12 hours before using.**